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
Fortgeschrittenen Praktikum

Teil I

Szintillationszähler

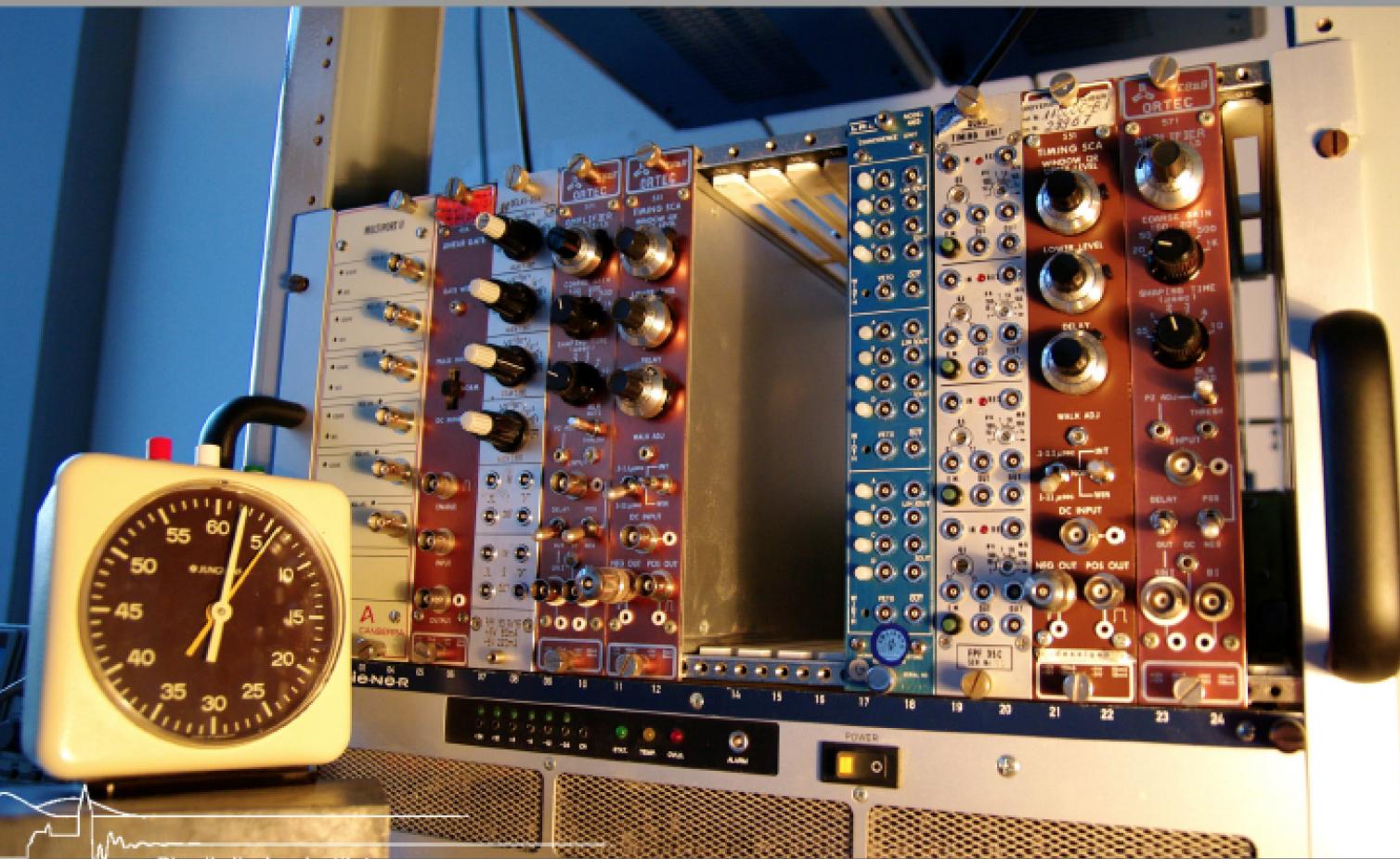FORTGESCHRITTENEN PRAKTIKUM I

TH. KLAUDOR-KLEINGROTHAUS (7/2016)

Szintillationszähler

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1 Introduction

In this experiment the decay spectra of different isotopes and the angle dependency of the pair annihilation will be measured. In parallel the student should get a deeper understanding in analog electronics, scintillators and NIM standards.

In the first part of the experiment the energy calibration of the readout system will be done with known radioactive sources and their spectra. For this the isotopes ^{22}Na , ^{152}Eu , ^{60}Co and ^{228}Th will be provided. The first three will be used for the energy calibration. Based on this calibration the spectrum of Thorium will be evaluated and compared to literature values.

In the second part the angle dependency of the pair annihilation will be measured by using the coincidence method. In addition to the NaI scintillator, which is used in the first part of the experiment, an organic scintillator will be used (??). The isotope ^{22}Na will be used, which decays via beta decay. For this part the different properties of the two kinds of scintillators have to be considered.

2 Preparation

Please ensure you're familiar with the provided devices, their functionality and with the provided isotopes. In the section ?? almost all information is provided which is necessary for this experiment. Additional information can be found in the datasheets of the manufacturers in chapter ??.

As in this experiment different NIM modules will be used ??, the basic functions and properties should be known at the beginning of the experiment. As preparation for the experiment the figures ?? und ?? have to be filled out. Figure ?? shows the circuit of the first part of the experiment and figure ?? the circuit of the second part. At each module the in- and outputs are indicated. Using the information from chapter ?? the corresponding signal shape of the signals can be put into the blank fields. Afterwards the modules have to be connected in a reasonable way for the experiment. In respect to the coincidence measurement there is a signal which has to be delayed.

In addition a general preparation to the topic is necessary. Questions might be:

- What are the characteristic signals of the provided isotopes?
- How do the niveau schemata look like? Have a look in the Staatsexamensarbeit or in chapter ??.
- Which decays are most probable and which peaks do you expect to see in the spectra?
- Which roles play α , β und γ decay in this experiment and how do they interact with the equipment?
- What does coincidence mean in this context?
- Think about other applications of the coincidence method in physics.

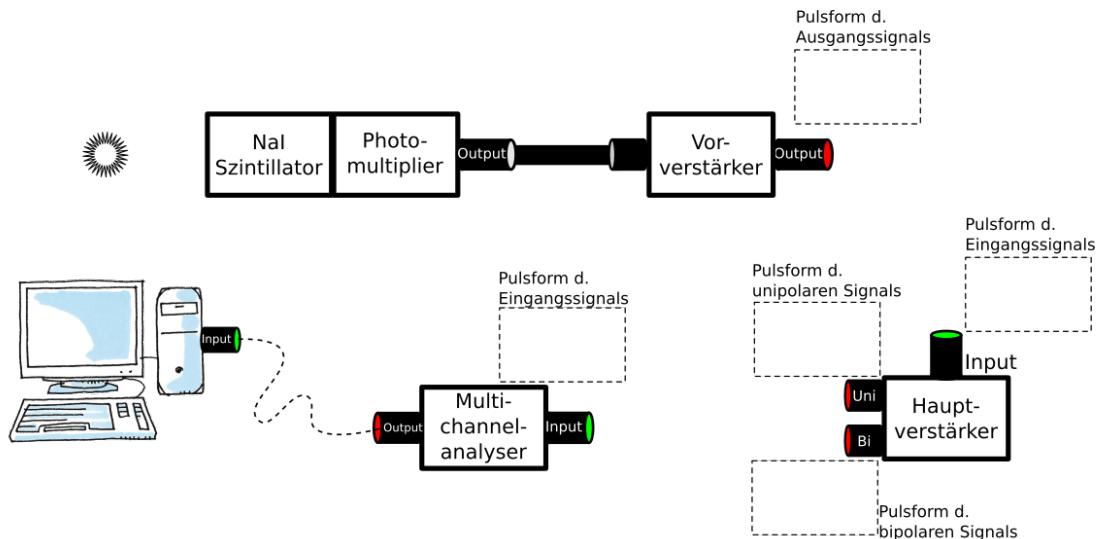
Abbildung 2: Circuit of part 1

3 Execution

Please check yourself to do all necessary measurements to answer the questions in chapter ??.

3.1 Part 1

1. Display the signals of the different outputs of the devices on the oscilloscope. Compare and verify them with the expected signals from figure ?? and ?? . Please consider chapter ?? — Signalabschluss) as it is important to get the experiment done.
2. Modify the settings of the devices and document their effect on the oscilloscope.
Place a radioactive probe at a reasonable position.
3. Modify the amplification in such way that the whole spectrum of the MCA is used for the measurements. As the isotopes emits different energies, the setting should be chosen such way that all peaks of each isotope can be seen. The amplification should not be modified after starting the measurements. Otherwise the energy calibration could not be done. How to readout the MCA with the PC can be looked up in section ??.
4. The measurement duration for the isotopes Europium, Cobalt and Natrium should be 30 min (minimum).

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5. The spectrum of Thorium should be measured over night.
6. The measurement of the background should take 3 hours.
7. Consider the comments for each part of the experiment in chapter ??.

3.2 Part 2

1. Set up the circuit as prepared in figure ??.
2. Discriminate the signal of the NaI scintillator for the corresponding energy of the process. Use your energy calibration! Put the radioactive at a reasonable position to ensure a high rate to make discrimination more easy.
3. Verify the amplification of the organic scintillator and correct it. The probe should be placed correspondingly. Why could the organic scintillator not be discriminated on a certain energy?
4. Place the probe to its position for the coincidence measurement.
5. Use the additional outputs of the timing unit to display the signals on the oscilloscope.
6. Use a reasonable display mode at the oscilloscope (section ??) to display the frequency of the occurrence of the signals to be observed.
7. Optimize the time delay of the ...scintillator to achieve a coincidence. In addition to the visual information of the oscilloscope, the HEX counter should be used to set an optimal time delay in a measurement series.
8. Decide on appropriate measurement times and angles in which you want to measure the angle dependence of pair annihilation.

4 Auswertung

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1. Prepare a circuit diagram.
2. Describe and discuss the observed signals (time, duration, pulse height, ...). Insert sketches or screenshots of the signals from the oscilloscope.

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571

Amplifier

Due to its versatility, the ORTEC Model 571 Amplifier is ideally suited for use with germanium detectors, silicon charged-particle detectors, proportional counters, scintillation counters, and pulsed ion chambers.

The 571 incorporates an automatic gated baseline restorer that offers exceptional performance over a very wide range of input counting rates (Fig. 1). The gated baseline restorer (BLR) includes a discriminator that operates the sensing circuits that normally establish the baseline reference for the MCA. Performance of the spectrometer depends on the precision of the setting of the BLR threshold. The Model 571 offers the convenience of an automatic threshold control, which typically gives as good or better results than those the most experienced operator could achieve manually.

The active filter networks permit the Model 571 to generate very symmetrical unipolar outputs with optimum signal-to-noise ratios over a wide range of time constants. The unipolar output can be internally delayed by 2 μ s via a front-panel control when desired for gating

purposes. The instrument also provides a bipolar output for timing and gating applications.

Any dc drift in an amplifier output causes spectrum broadening. The excellent dc stability of the Model 571 eliminates spectrum broadening caused by dc drift and ensures that the high resolution capability of germanium detectors is realized.

Fig. 2. BLR Auto Control.

- General-purpose amplifier for energy spectroscopy with all types of detectors
- Unipolar and bipolar outputs
- Selectable unipolar output delay
- Active filter networks with wide range of time constants
- Wide gain range
- Gated baseline restorer with automatic BLR threshold control for excellent counting rate performance

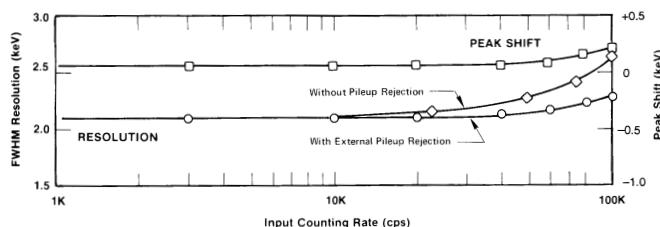
Fig. 1. Typical Resolution and Baseline Stability vs. Counting Rate for the Model 571 in a γ -Spectroscopy System.

3.29

571

Amplifier

Specifications

PERFORMANCE

GAIN RANGE Continuously adjustable from 1 to 1500.

PULSE SHAPE Semi-Gaussian on all ranges with peaking time equal to 2.2τ and pulse width at 0.1% level equal to 2.9 times the peaking time.

INTEGRAL NONLINEARITY For 2- μ s shaping time, $<\!0.05\%$.

NOISE Typically $<5 \mu\text{V}$ for unipolar output referred to the input, using 2- μ s shaping and Coarse Gain ≥ 100 .

TEMPERATURE INSTABILITY

Gain $\leq 0.0075\%/\text{^{\circ}C}$, 0 to 50 $^{\circ}\text{C}$.

DC Level $\leq 50 \mu\text{V}/\text{^{\circ}C}$, 0 to 50 $^{\circ}\text{C}$.

BIPOLAR CROSSOVER WALK $\leq \pm 3$ ns at 0.5 μ s for 50:1 dynamic range, including contribution of an ORTEC Model 552 Single-Channel Analyzer.

OVERLOAD RECOVERY Recovers to within 2% of rated output from X800 overload in 2.5 nonoverloaded pulse widths using maximum gain for unipolar output. Same recovery from X1000 overload for bipolar.

SPECTRUM BROADENING Typically $<16\%$ FWHM for a ^{60}Co 1.33-MeV gamma line at 85% of full scale for an incoming count rate of 1 to 100,000 counts/s. Unipolar output, 2- μ s shaping.

SPECTRUM SHIFT Peak position shifts typically $<0.024\%$ for a ^{60}Co 1.33-MeV gamma line at 85% of full scale, measured from 1 to 100,000 counts/s at the unipolar output, 2- μ s shaping.

CONTROLS

FINE GAIN 10-turn precision potentiometer with graduated dial for continuously variable direct-reading gain factor X0.5 to X1.5.

COARSE GAIN 6-position switch selects feedback resistors for gain factors of 20, 50, 100, 200, 500, and 1k. Jumper on the printed wiring board (PWB) selects X0.1 attenuation.

SHAPING TIME 6-position switch selects time constants for active pulse-shaping filter network from 0.5, 1, 2, 3, 6, and 10 μ s.

DELAY Locking toggle switch selects either 2- μ s delay (In) or prompt (Out) timing for unipolar output.

INPUT POLARITY Locking toggle switch selects either Pos or Neg input pulse polarity.

PZ ADJ Screwdriver adjustable potentiometer to set the pole-zero cancellation to compensate input decay times from 40 μ s to ∞ .

BLR 3-position locking toggle switch selects the source of control for the gated baseline restorer discriminator threshold from:

Auto The BLR threshold is automatically set to an optimum level, as a function of the signal noise, by an internal circuit.

PZ Adj The BLR threshold is determined by the threshold potentiometer. The BLR time constant is also greatly increased to facilitate PZ adjustment; this position may give the lowest noise for count rates under 5000 counts/s and/or longer shaping times.

Threshold The BLR threshold is manually set by the threshold potentiometer.


DC Screwdriver adjustable potentiometer to set the unipolar output dc level; range ± 100 mV.

INPUT


INPUT BNC front- and rear-panel connectors accept either positive or negative pulses with rise time of 10 to 650 ns and decay times of 40 μ s to ∞ ; $Z_{in} \leq 1000 \Omega$ dc-coupled; linear maximum 10 V; absolute maximum 20 V.

OUTPUTS

UNI Front-panel BNC connector with $Z_o < 1 \Omega$ and rear-panel connector with $Z_o = 93 \Omega$, short circuit proof; prompt or delayed with full-

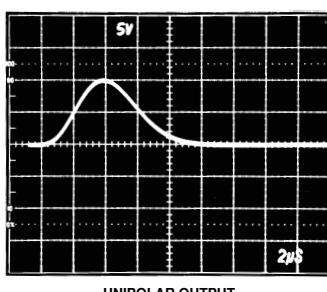UNIPOLAR OUTPUT

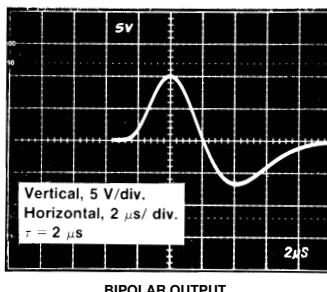BIPOLAR OUTPUT

scale linear range of 0 to +10 V; active filter shaped; dc-restored, dc level adjustable to ± 100 mV.

BI Front-panel BNC connector with $Z_o < 1 \Omega$ and rear-panel connector with $Z_o = 93 \Omega$, short circuit proof; prompt output with positive lobe leading and linear range of ± 10 V; active filter shaped.

PREAMP POWER Rear-panel standard ORTEC power connector, Amphenol 17-10090, mates with captive and noncaptive power cords on all ORTEC preamplifiers.

BUSY OUTPUT Rear-panel BNC connector with $Z_o < 10 \Omega$ provides a +5 V logic pulse for the duration that the input pulse exceeds the baseline restorer discriminator.

ELECTRICAL AND MECHANICAL

POWER REQUIRED +12 V, 60 mA; -12 V, 30 mA; +24 V, 100 mA; -24 V, 105 mA.

WEIGHT

Net 1.5 kg (3.3 lb).

Shipping 3.1 kg (7.0 lb).

DIMENSIONS Standard single-width NIM module 3.43 X 22.13 cm (1.35 X 8.714 in.) per DOE/ER-0457T.

Ordering Information

To order, specify:

Model Description

571 Amplifier

7.2 Single Channel Analyser

Single-Channel Pulse-Height Analyzers

Contents

Introduction	4.1
Single-Channel Analyzer Applications Guide	4.3
Single-Channel Analyzer Selection Guide	4.4
550A Single-Channel Analyzer	4.5
551 Timing Single-Channel Analyzer	4.7
552 Pulse-Shape Analyzer/Timing SCA	4.9
583A Constant-Fraction Differential Discriminator/SCA	6.16
590A Amplifier and Timing Single-Channel Analyzer	3.40
850 Quad Single-Channel Analyzer	4.12

Introduction

The amplitude of the analog pulse at the output of a spectroscopy amplifier is typically proportional to the charge released in the detector or to the energy of the detected event. Selection of a range of signal levels at the output of the amplifier is equivalent to the selection of a range of energies or charge for these events. This selection can be accomplished by the use of discriminators and single-channel analyzers (SCAs). A discriminator produces an output logic pulse only if its input signal exceeds a preset threshold level. A single-channel analyzer produces an output logic pulse only if the peak amplitude of its input signal falls within the pulse-height window that is established with two preset threshold levels.

Figure 1 shows three pulses that might be provided from a main amplifier to an integral discriminator. The first pulse has an amplitude less than the adjusted discriminator threshold and generates no output logic signal. Each of the last two pulses has sufficient amplitude to produce an output logic signal. The output signals indicated in Fig. 1 are generated when the leading edge of the input signal crosses the

discriminator threshold level. Therefore, the time of the output response is a function of the amplitude and rise time of the input signals. This amplitude and rise time dependence leads to "time walk" of the output signal relative to the beginning of the input pulse. The discriminator output is produced earlier by pulses with larger amplitudes and later by pulses with lower amplitudes.

Figure 2 shows three pulses that might be provided from a main amplifier to an SCA. Only the B pulse satisfies the conditions necessary to produce an SCA output logic signal.

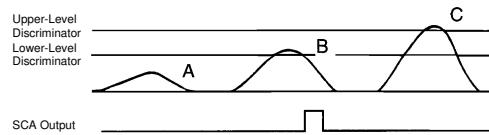
Fig. 2. Single-Channel Analyzer Function.

Removal of the upper-level-discriminator restrictions from the SCA allows it to be used as an integral discriminator. If the upper-level restrictions were removed from the unit whose output is shown in Fig. 2, both pulses B and C would be marked by logic outputs.

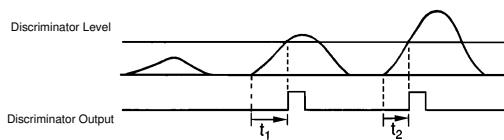
Fig. 1. Integral Discriminator Output Triggering.

Single-Channel Pulse-Height Analyzers

Three primary modes of discriminator operation are available in ORTEC SCAs: Integral, Normal, and Window. In the Integral mode of operation, the SCA can function as an integral discriminator, as indicated in the preceding paragraph. In the SCA Normal mode of operation, the upper-level and lower-level thresholds are independently adjustable. In the SCA Window mode, the upper-level threshold control is used to establish a voltage level that is added to the lower-level threshold voltage to yield the upper-level discriminator (ULD) threshold level. Thus, when the lower-level setting is changed, the upper-level threshold changes by the same amount. An external voltage reference for the lower-level discriminator (LLD) can be supplied to scan the window through a preselected range of pulse heights.

Unlike an integral discriminator, the output logic signal from a single-channel analyzer must be produced after the input pulse reaches its maximum amplitude. This timing sequence must provide sufficient time for the SCA logic circuitry to determine if the input signal exceeded the upper-level threshold.

ORTEC provides two basic types or classifications of SCAs: non-timing SCAs and timing SCAs. The technique used to produce the output logic signals from an SCA determines its classification. Non-timing units, such as the Models 550A, and 850, produce an SCA output pulse if the input signal is within the window settings. The output occurs when the trailing edge of the input signal recrosses the lower-level threshold. Figure 3 shows two superimposed output pulses from a main amplifier that meet the window requirements of the single-channel analyzer. The output from the non-timing SCA for each pulse is shown below the pulses. Since the linear input pulses are referenced to the same starting time, it is clear that the output logic signals exhibit "time walk" relative to the input pulses.

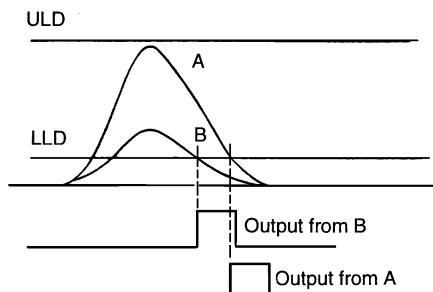
Fig. 3. Non-Timing SCA Output Triggering.

Timing SCAs, such as the ORTEC Models 551, 552, and 590A, produce SCA output logic signals that are precisely related in time to the occurrence of the event being measured. This time relationship implies that the time of occurrence of the SCA output signal is "walk-free" or nominally independent of the amplitude of the input signal, for a given rise time. In addition to simple counting applications, the time-related output can be used for coincidence measurement, pulse-shape discrimination, and other applications where the precise time of occurrence is important.

Figure 4 shows two pulses from a main amplifier and the response for a peak-detection single-channel analyzer such as the Model 590A Amplifier and Timing Single-Channel Analyzer. Although the amplitudes of the amplifier pulses differ, their peaks occur at approximately the same time, and the SCA outputs are produced when the peaks of the input pulses are detected.

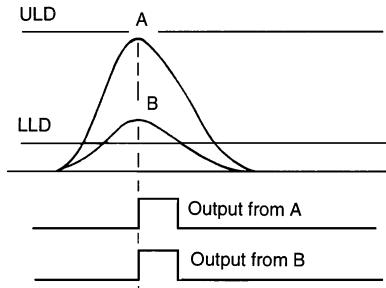
Fig. 4. Peak-Sensing SCA Output Triggering.

The conventional zero-crossing technique has been widely used for timing single-channel analyzers. This technique utilizes the zero-crossing of the bipolar output signal from a pulse-shaping amplifier to derive timing information, and uses the peak amplitude of the pulse for the energy range information. Figure 5 shows two bipolar pulses provided from a main shaping amplifier. Both pulses meet the SCA window requirements. Each output signal is generated when the corresponding input signal crosses the baseline. Figure 5 illustrates that the time of occurrence of the SCA output signals is precisely related to the occurrence of the detected event and is independent of input signal amplitude. Either double-delay-line-shaped pulses or RC-shaped pulses may be used, but the former provide better timing resolution. The bipolar output from delay-line amplifiers such as the Model 460 is well suited to zero-crossover timing with the ORTEC Model 552, because the input signal crosses the baseline with a large slope even when the pulse amplitude is low.

Single-Channel Pulse-Height Analyzers

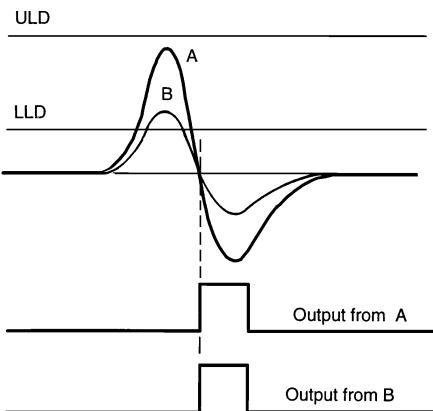
Fig. 5. Zero-Crossover SCA Output Triggering.

The bipolar output signal from a double-delay-line shaping amplifier crosses the baseline at a fixed fraction that is effectively 50% of the charge collected from the detector. Thus, conventional zero-crossing timing can be considered as timing at a constant fraction of the input signal amplitude. A trailing-edge constant-fraction technique* can be used with either unipolar or bipolar signals to derive a time-pickoff pulse after the peak time of the signal from the shaping amplifier. This technique is extremely useful when incorporated in timing single-channel analyzers. Figure 6 illustrates the trailing-edge constant-fraction technique for two unipolar input signals of identical rise times but different amplitudes. The time of occurrence of the output signals is independent of output signal amplitudes.

The trailing-edge constant-fraction timing technique is available with two ORTEC SCAs: Models 551 and 552.

The Model 552 can also be used as a pulse-shape analyzer. The best known application of this technique is in the separation of the neutron and gamma responses of some scintillators (see page A.18). Collection time differences for the two types of radiation result in shape or rise time variations in the signals from a spectroscopy amplifier. When used with an ORTEC Time-to-Amplitude Converter, the Model 552 can resolve these shape variations over a 200:1 dynamic range of input signal amplitudes. The Model 552 accomplishes the shape measurement of the input signals by evaluating the timing at two different fractions.

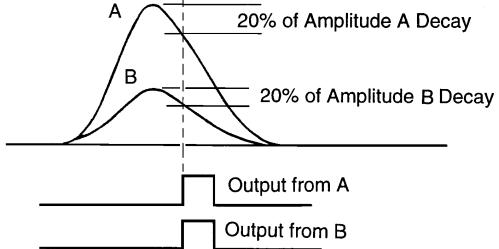
Fig. 6. Constant-Fraction SCA Output Triggering.

The SCA function can also be applied to fast analog signals in rising-edge constant-fraction discriminators. See the Fast-Timing Discriminators Section for details.

The following Selection Guide provides comparative data for all ORTEC Single-Channel Analyzers.

Single-Channel Analyzer Applications Guide

Model	Recommendations
550A	Versatile, economical, general-purpose counting.
551	SCA plus constant-fraction timing.
552	SCA plus constant-fraction timing and pulse-shape analysis.
583A	SCA function in a fast-timing constant-fraction discriminator for signals from a photomultiplier anode.
590A	Cost efficient, includes built-in amplifier.
850	Economical, four SCAs in a single-width module for general-purpose counting.

*The basic circuit for implementing this technique is patented by ORTEC, U.S. Patent No. 3,714,464.

Timing Single-Channel Analyzer

The ORTEC Model 551 Timing Single-Channel Analyzer performs the dual functions of single-channel pulse-height analysis and timing signal derivation.

The patented* trailing-edge constant-fraction timing technique provides unexcelled timing on either unipolar or bipolar signals and shows better results than are possible with conventional leading-edge discriminators.

With SCAs that utilize leading-edge timing, the rise time of the input pulses causes degradation of time resolution because the pulses have varying amplitudes.

Constant-fraction timing compensates for varying amplitudes and essentially eliminates this timing shift, giving consistently better timing results.

For the internally set 50% fraction, the output occurs soon after the midpoint on the linear input trailing edge to facilitate gating and accumulation of data at very high input rates. This technique also minimizes timing shift and dead time when used with sodium iodide, silicon, and germanium detectors, thereby allowing better system time resolution and higher counting rates.

The constant-fraction technique makes it possible to realize significant improvements in time resolution in most timing applications. Notice that analysis is made of the main amplifier output. This technique allows optimization of time resolution and extension of dynamic range for neutron-gamma discrimination and other timing applications. Walk of <3 ns for 100:1 dynamic range using input pulses from a pulser is possible.

The Model 551 is versatile, with three basic operating modes provided. In the Window mode, the unit operates as a high-resolution, narrow (0 to 10%) window, single-channel analyzer. For wide-window applications, the Normal mode is used. In this mode the upper-level and lower-level controls are independently variable from 0 to 10 V, and an output is generated for pulses analyzed between the levels. Through use of the separate rear-panel LL Out

- Single-channel analyzer and timing signal derivation
- Trailing-edge constant-fraction timing provides walk $<\pm 3$ ns for 100:1 dynamic range
- Integral, normal, and window modes
- Separate lower-level and upper-level discriminator outputs
- DC-coupled
- Adjustable delay 0.1 to 11 μ s
- Provision for external baseline sweep

and UL Out outputs, the unit can operate as a dual wide-dynamic-range integral discriminator for leading-edge timing or for pulse routing.

The dc-coupled input of the Model 551 makes it possible to take full advantage of the baseline restoration of the main amplifier for maximum performance at widely varying counting rates.

The continuously adjustable output delay (two ranges covering 0.1 to 11 μ s) makes it possible to align output signals that have actual time differences without a need for additional delay devices or modules. Alternatively an External strobe input can be used to cause an SCA output at the desired time.

For an application where it is desirable to scan an entire spectrum, an external base-line sweep input is provided via the rear-panel LL Ref Ext BNC connector. In this mode of operation, the baseline (lower-level threshold) on which a window is riding is swept through an energy range and the count rate is recorded as a function of energy.

*U.S. Patent No. 3,714,464.

Timing Single-Channel Analyzer

Specifications

PERFORMANCE

DYNAMIC RANGE 200:1.

PULSE-PAIR RESOLVING TIME Output pulse width plus Delay (as selected by the front-panel Delay controls), plus 100 ns for fast NIM output or plus 200 ns for positive NIM output. Minimum resolving time for negative output 220 ns; for positive output 800 ns.

THRESHOLD TEMPERATURE

INSTABILITY $\leq \pm 0.01\%/\text{C}$ of full scale, 0 to 50°C using a NIM Class A power supply (referenced to -12 V).

DISCRIMINATOR NONLINEARITY $\leq \pm 0.25\%$ of full scale (integral) for both discriminators.

DELAY TEMPERATURE INSTABILITY $\leq \pm 0.03\%/\text{C}$ of full scale, 0 to 50°C.

DELAY NONLINEARITY $< \pm 2\%$ of delay range.

WINDOW WIDTH CONSTANCY $\leq \pm 0.1\%$ variation of full-scale window width over the linear range 0 to 10 V.

MINIMUM INPUT THRESHOLD 50 mV for lower-level discriminator.

TIME SHIFT vs PULSE HEIGHT (WALK)

Walk (ns)		Dynamic Range
System A	System B	
±1.0	±2.0	10:1
±2.5	±4.0	50:1
±3.0	±4.0	100:1

System A: Using an ORTEC Model 460 Amplifier, single delay-line mode, integrate $\leq 0.1 \mu\text{s}$ with 1- μs delay line.

System B: Using an ORTEC Model 570, 571, or 572 Amplifier, unipolar output with 0.5- μs shaping time. Input from ORTEC Model 419 Pulser.

CONTROLS

LOWER LEVEL Front-panel 10-turn potentiometer adjustable from 0 to 10 V; when the rear-panel LL Ref mode switch is set on Int, determines the threshold setting for the lower-level discriminator. When the LL REF mode switch on the rear panel is in the EXT position, this control is ineffective.

WINDOW OR UPPER LEVEL Front-panel 10-turn potentiometer determines the window width (0 to +1 V) in the Window mode or the upper-level (0 to +10 V) threshold in the Normal mode. This control is disabled in the Integral mode.

INT/NOR/WIN Front-panel 3-position locking toggle switch selects one of three operating modes:

Integral LL sets a single-discriminator threshold (0 to +10 V) and UL is disabled.

Normal UL and LL are independently adjustable levels (0 to +10 V).

Window LL sets the baseline level (0 to +10 V) and UL sets the window width (0 to +1 V).

DELAY RANGE Front-panel locking toggle switch selects delay ranges of 0.1 to 1.1 μs or 1.0 to 11 μs .

DELAY Front-panel 10-turn potentiometer for continuous adjustment of output delay over selected range. In the external strobe mode the delay control adjusts the automatic reset time from $\approx 5 \mu\text{s}$ to 50 μs .

WALK ADJUST Front-panel screwdriver adjustment for precise setting of walk compensation.

LL REF MODE Rear-panel 2-position locking toggle switch selects either the front-panel LL potentiometer or the voltage signal applied to the rear-panel LL REF EXT connector as the LL discriminator reference threshold.

STROBE Rear-panel 2-position locking toggle switch selects either Internal or External source for the SCA output signal strobe function.

INPUTS

SIGNAL INPUT Front-panel dc-coupled BNC connector accepts positive unipolar or bipolar signal, 0 to +10 V linear range, $\pm 12 \text{ V}$ maximum; width 100 ns; $1000\text{-}\Omega$ input impedance. Rear-panel ac-coupled BNC connector accepts positive unipolar or bipolar signal, 0 to +10 V linear range, $\pm 100 \text{ V}$ maximum; width 0.2 to 10 μs ; $1000\text{-}\Omega$ input impedance.

LL REF EXT When the rear-panel LL REF mode switch is on EXT, the rear-panel LL REF EXT BNC connector accepts the lower-level biasing (an input of 0 to -10 V on this connector corresponds to a range of 0 to 10 V for the lower-level discriminator setting). Input protected to $\pm 24 \text{ V}$.

EXT STROBE INT When the rear-panel EXT/INT STROBE locking toggle switch is in EXT, the rear-panel EXT STROBE IN BNC connector accepts a positive NIM-standard input, nominally +5 V, 500 ns wide, to cause an output to occur from the SCA. The external strobe should be given within 5 μs (or 50 μs as determined by the front-panel Delay control) of the linear input. At the end of this period, the Model 551 resets its internal logic without producing an output signal.

OUTPUTS

SCA POS OUT Front- and rear-panel BNC connectors provide positive NIM-standard output, nominally +5 V; 500 ns wide; $10\text{-}\Omega$ output impedance. For internal strobe the output occurs at the midpoint of the linear input trailing edge plus the output Delay as selected by the front-panel controls. For external strobe the output occurs at the time of strobe signal.

SCA NEG OUT Front-panel BNC connector provides fast NIM-standard output, nominally -16 mA (-800 mV on $50\text{-}\Omega$ load); width $\leq 20 \text{ ns}$; rise time $\leq 5 \text{ ns}$; $\leq 10\text{-}\Omega$ output impedance. Output occurs at the mid-point of the linear trailing edge plus the output Delay as selected by the front-panel controls.

LL OUT Rear-panel BNC connector provides positive NIM-standard output, nominally +5 V, 500 ns wide; $\leq 10\text{-}\Omega$ output impedance. Output occurs as leading edge of linear input crosses the LL threshold.

UL OUT Rear-panel BNC connector provides NIM-standard output, nominally +5 V, 500 ns wide; $\leq 10\text{-}\Omega$ output impedance. Output occurs as leading edge of linear input crosses the UL threshold.

ELECTRICAL AND MECHANICAL

POWER REQUIRED +12 V, 160 mA; -12 V, 110 mA; +24 V, 90 mA; -24 V, 50 mA.

WEIGHT

Net 1.1 kg (2.5 lb).

Shipping 2.25 kg (5.0 lb).

DIMENSIONS NIM-standard single-width module 3.43 X 22.13 cm (1.35 X 4.714 in.) per DOE/ER-0457T.

Related Equipment

The Model 551 is compatible with all ORTEC amplifiers and other amplifiers having a 0 to 10 V positive, linear output range.

Ordering Information

To order, specify:

Model	Description
551	Timing Single-Channel Analyzer

7.3 Multiport

Features

- Up to six independent Multichannel Analyzers (in-the-field customer upgradable)
- Each MCA fully computer controlled
- 16K ADC, 1 μ s fast fixed conversion time with linearization enhancement circuit for excellent non-linearity
- Fully buffered memory allowing simultaneous ADC memory write and spectrum transfer to host computer
- PHA and MCS acquisition modes
- Full I/O control (including sample changer)
- On-board data backup
- Suitable for HPGe, NaI(Tl), CdTe, PIPS®, X-PIPS™ and other detector technologies
- Communications interface: Ethernet/USB or USB only
- Fully supported by Genie™ 2000 software

Multiport II Multichannel Analyzer

Description

The Multiport II is a double width NIM intended for use where existing amplifier and HVPS modules are to be combined with a high analog performance, low-cost ADC/MCA. It provides greater flexibility than traditional bus plug-in boards but at a comparable cost. Multiport II can be purchased in versions from one to six inputs. Units with fewer than six inputs can be field upgraded with additional inputs. Two communications interface versions are available: Ethernet/USB and USB only.

The Multiport II is easy to install and easy to link. Host computer interfacing is accomplished via a standard USB port or a standard Ethernet port (unshielded twisted pair) and through standard protocols. These widely recognized interfaces make the Multiport II compatible with a wide range of computer platforms, so the Multiport II user need not worry about compatibility as computer bus standards evolve. Also, multiple Multiport II units can be totally remote-controlled from a single computer. The Multiport II provides better performance and lower noise than plug-in PC board based MCAs by keeping all sensitive components out of the computer itself and by putting the MCA close to the radiation detectors. Furthermore, there is only one cable between the host computer and the Multiport II NIM module for simplicity of inter connection. Multiport II is suitable for use with a wide range of radiation detectors. Selecting a proper preamplifier, amplifier and high voltage power supply, makes Multiport II compatible with NaI(Tl), HPGe, SiLi, CdTe, ion implanted, plastic scintillation, BGO and other detector technologies. Modular NIM packaging makes it easy to reconfigure systems as needs change or as new technologies become available.

The use of parallel processing technique results in an effective zero time to "add one" in memory, resulting in a fast conversion time. Furthermore, the on-board implementation of the sliding scale method improves the linearity and the channel uniformity. To protect long collections from power failure, the Multiport II features on-board data and settings retention at fixed time intervals for each of its MCAs.

Multiport II provides full I/O support, including Canberra standard PUR/LTC, sample changer synchronization and advanced PHA operations.

For more flexibility and more applications, Multiport II supports both Pulse Height Analysis (PHA) and Multichannel Scaling (MCS) modes of operation. In MCA mode, a single channel analyzer (SCA) is also available as well as input/output that allows acquisition to be synchronized by external modules.

Multiport II Multichannel Analyzer

SCA output works in both modes, PHA and MCS. In MCS mode, SCA pulses may also be selected as input pulses. The maximum input rate is 10 MHz in TTL input mode.

Software Support

The Multiport II is fully supported by Genie 2000 software family suite. Multiport II is fully remote controlled under Genie 2000 (parameters such as gain, LLD, ULD, ADC zero and MCS parameters) via standard Genie 2000 tools. Genie 2000 supports a wide range of time proven spectral analysis algorithms, modern spectrum display and user interface as well as a variety of special applications programs. Genie 2000 software solutions are available for applications ranging as widely as laboratory gamma and alpha spectroscopy, waste measurement, whole body counting and nuclear safeguards.

Specifications

ADC INPUT

- ADC IN – Accepts positive unipolar or bipolar, positive lobe leading, pulses for PHA and MCS (Internal Pulse Source); amplitude 0 to +10 V or 0 to +5 V; switch selectable; rise time >0.25 μ s; width >0.5 μ s; Z_{in} = 1.33 k Ω (10 V range) or 2 k Ω (5 V range); direct coupled; front and rear panel BNCs, one pair per MCA.

LOGIC INPUTS

Rear panel mini-D I/O connectors for logic input and output signals, one per MCA.

- XINPUT – External start/stop and sample changer ready input; active high, active low setting; 10 k Ω pull-up to 5 V; TTL pulse; pin 1 of the I/O connector.
- XSUSPEND – Suspends pulse processing and preset counters; fixed active low; 10 k Ω pull-up to +5 V; TTL pulse; pin 3 of the I/O connector.
- XDT – PUR dead time input; active high, active low setting; forced to inactive state by 1 k Ω pull-up to +5 V or pull-down to GND when no signal is present; TTL pulse; pin 5 of the I/O connector.
- XMCSRESTART – MCS sweep advance input; active high, active low setting; edge triggered; 10 k Ω pull-up to +5 V; TTL pulse; pin 10 of the I/O connector.
- XMCSADVANCE – External channel advance input; active high, active low setting; edge triggered; duration \geq 10 ns; max. rate 1 MHz; 10 k Ω pull-up to +5 V; positive TTL pulse; pin 11 of the I/O connector.
- XMCSPULSE – MCS input; active high, active low setting; edge triggered; duration \geq 10 ns; rate \leq 10 MHz; 10 k Ω pull-up to +5 V; positive TTL pulse; pin 12 of the I/O connector.
- XCOINC – Coincidence/Anticoincidence input; active high, active low setting; edge triggered, level sensitive selection; forced to inactive state by 1 k Ω pull-up to +5 V or pull-down to GND when no signal is present; TTL pulse; pin 13 of the I/O connector.
- XREJ – Pileup reject input; must occur during the ADC linear gate (LG) time active high, active low setting; forced to inactive state by 1 k Ω pull-up to +5 V or pull-down to GND when no signal is present; positive TTL pulse; pin 14 of the I/O connector.

LOGIC OUTPUTS

Rear panel mini-D I/O connectors for logic input and output signals, one per MCA.

- XOUTPUT – Advance sample changer output; active high, active low setting; 150 ms TTL pulse; pin 2 of the I/O connector.
- XSCA – Single channel analyzer output; \approx 250 ns TTL pulse generated for input pulses detected between LLD and ULD; active high; pin 4 of the I/O connector.
- XLG – PUR linear gate output; active when ADC acquires an input pulse; active high, active low setting; TTL pulse; pin 6 of the I/O connector.
- XCOLLECTSTATUS – External acquire status; active high, active low setting; TTL pulse; pin 9 of the I/O connector.
- +5 V – Power output for external circuitry; 100 mA max. per connector; pin 7 of the I/O connector.

REAR PANEL CONTROL

- RESET – Resets the CPU and communication; stops any acquisition in progress; spectra and MCA status stored in flash memory not affected; rear panel pushbutton.

INTERNAL CONTROLS

- ADC INPUT RANGE – Sets the ADC input range to 0–10 V (default) or 0–5 V; slide switch on the MCA board.
- POWER SUPPLY VOLTAGE – Selects either \pm 12 V (default) or \pm 6 V; slide switch on the main board.

SOFTWARE CONTROLS

- ADC LLD – 0 to 100% of full scale (4096 steps).
- ADC ULD – 0 to 110% of full scale (4096 steps).
- ADC ZERO – \pm 2.5% of full scale (4096 steps).
- ADC CONVERSION GAIN – 256, 512, 1024, 2048, 4096, 8192, 16 384.
- REAL TIME PRESET – Maximum of \approx 49.71 days.
- LIVE TIME PRESET – Maximum of \approx 49.71 days.

INDICATORS

- MCA #n – ON when MCA board in slot #n is installed and power in ON; front panel green LED.
- ACQUIRE – ON when the corresponding MCA is acquiring; front panel green LED.
- RATE – Flashes for every input pulse processed in the corresponding MCA; front panel yellow LED.
- COMM – ON when unit is ready for USB and/or Ethernet communication; flashes brighter when data transfer occurs; front panel yellow LED.

COMMUNICATION

- USB – Standard USB port; rear panel USB type-B connector.
- ETHERNET – Optional Ethernet port; rear panel RJ-45 connector.

PERFORMANCE

- INTEGRAL NONLINEARITY – $<\pm$ 0.025% of full scale over the top 99.5% of selected range.
- DIFFERENTIAL NONLINEARITY – $<\pm$ 0.9% of full scale over the top 99.5% of range including effects from differential nonlinearity.
- GAIN DRIFT – $<\pm$ 0.005% of full scale/ $^{\circ}$ C.
- ZERO DRIFT – $<\pm$ 0.005% of full scale/ $^{\circ}$ C.

Multiport II Multichannel Analyzer

- LONG TERM DRIFT – <±0.005% of full scale/24 h at a constant temperature.
- PEAK SHIFT – <±0.025% of full scale at rates up to 100 kHz.
- ADC DEAD TIME – Linear gate time +1.2 µs.
- CHANNEL PROFILE – Typically flat over 90% of channel width.

MCS –

- INPUT RATE – Up to 10 MHz.
- DWELL TIME – Software selectable from 1 µs to ≈71 min, with varying resolution in increments of 1 µs.
- SWEEP SELECTION – Software selectable preset sweeps from 1 to 2^{32} –1.
- DEAD TIME – Between channels: 0. Between sweeps: 0 for internal sweep restart, 0–40 ns for external sweep restart.
- NUMBER OF CHANNELS – 16 384.
- COUNTS/CHANNEL – 2^{32} –1.

POWER

Depends on the setting of the internal control for Power Supply Voltage (± 12 V or ± 6 V):

±12 V Power Supply –

One MCA board:
+12 V dc – 250 mA –12 V dc – 50 mA

Each additional MCA board:
+12 V dc – 90 mA –12 V dc – 50 mA

±6 V Power Supply –

One MCA board:
+6 V dc – 475 mA –6 V dc – 50 mA

Each additional MCA board:
+6 V dc – 150 mA –6 V dc – 50 mA

PHYSICAL

- SIZE – Standard double width NIM module 6.86 x 22.12 cm (2.70 x 8.71 in.) per DOE/ER-0457T.
- NET WEIGHT – 1.36 kg (3.0 lb) with one MCA/ADC board.

ENVIRONMENTAL


- OPERATING TEMPERATURE – 0–45 °C.
- OPERATING HUMIDITY – 0–80% relative, non-condensing. Meets the environmental conditions specified by EN 61010, Installation Category I, Pollution Degree 2.

SYSTEM REQUIREMENTS

- A PC-based system with a USB port. Ethernet version also requires Ethernet connection.
- Model S500 Genie 2000 Basic Spectroscopy Software, V2.1a or later.

Ordering Information

- MP2- n U – Multiport II USB module with n inputs, where n is 1 to 6; not field upgradable to the Ethernet module; includes one 3 m shielded USB cable.
- MP2- n E – Multiport II Ethernet/USB module with n inputs, where n is 1 to 6; includes one 3 m shielded USB cable.
- MPT2-MCA – One additional Multiport II MCA/ADC board.
- C1801 – Multiport II PUR/LTC Cable (included with each MCA/ADC board).
- C1802 – Multiport II PUR/LTC Start/Stop and Sample Changer Cable.
- C1804 – Multiport II PUR/LTC and MCS Start/Stop Cable.

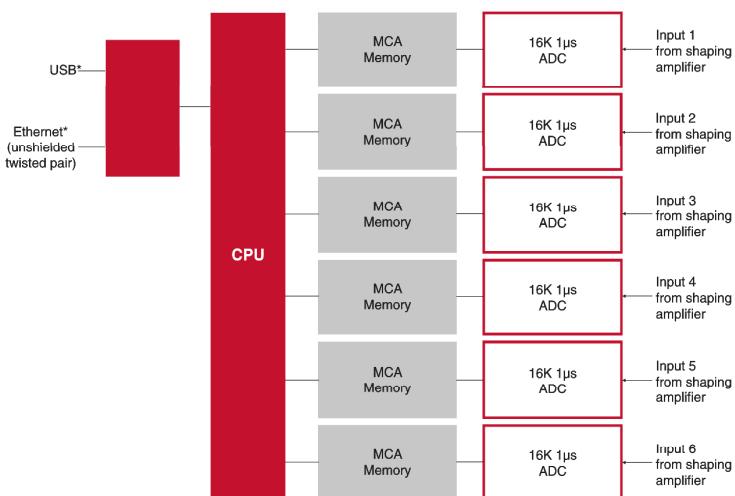*Depending on communication interface version.


X-PIPS and Genie are trademarks of Canberra Industries, Inc.
PIPS is a registered trademark of Canberra Industries, Inc.

©2008 Canberra Industries, Inc. All rights reserved.

7.4 Linear Gate

426 Linear Gate

- For passing and blocking analog signals in the range from +0.2 to +10 V
- Ungated or gated with coincidence or anticoincidence gating
- External or internal control of gate pulse width

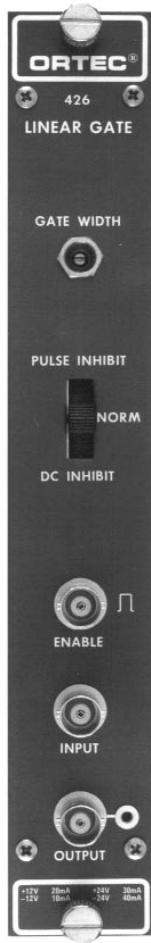The ORTEC Model 426 Linear Gate provides a variable gate duration with width controlled by a single-turn front-panel-mounted potentiometer. The nominal gate duration is from 0.3 to 4 μ s. Operation of the linear gate is controlled by a positive enable pulse. It is useful for selecting or inhibiting linear signals according to chosen coincidence or timing requirements.

The ORTEC Model 426 has two operating modes: all input signals not accompanied by an enable pulse are blocked or all signals are passed unless accompanied by an inhibit signal. The inhibit signal can be fed into the front-panel Enable connector for Pulse Inhibit operation or into the DC Inhibit connector for dc or continuous inhibit operation. The DC Inhibit mode provides external control of the gating period.

Specifications

PERFORMANCE

GAIN

Unity.
INTEGRAL NONLINEARITY $<\pm 0.15\%$ from 0.2 to 10 V.

PULSE FEEDTHROUGH <10 mV with a 10-V input pulse.

TEMPERATURE INSTABILITY $<\pm 0.015\%/\text{C}$, 0 to 50 $^{\circ}\text{C}$.

COUNTING RATE The gain shift of a 4-V reference pulse is $<0.25\%$ with the application of an additional count rate of 65,000 counts/s of 6-V random pulses.

CONTROLS

GATE WIDTH Continuously variable from 0.3 to 4 μ s.

OUTPUT PEDESTAL Adjustable to <1 mV.

PULSE INHIBIT/NORM/DC INHIBIT 3-position mode switch permits selection of the function of any pulse or dc level furnished through the front-panel Enable Input connector, or the rear-panel DC Inhibit connector.

Norm Input pulse will be gated through to the output during a gate width interval following the leading edge of each Enable Input pulse.

Pulse Inhibit Input pulses will be inhibited from passing through the output during a gate width interval following each Enable Input pulse.

DC Inhibit Input pulses will be inhibited from passing through the output during intervals of pulses or dc levels through the rear-panel DC Inhibit connector.

INPUTS

LINEAR INPUT Unipolar or bipolar with positive portion leading. Rated range 0.2 to 10 V, 12 V maximum. Input impedance >5000 Ω ; BNC connector. Input is ac-coupled with a passive symmetric baseline restorer. BLR can be bypassed for dc-coupling.

ENABLE OR INHIBIT INPUT Any positive input >2 V, maximum input 20 V. Enable impedance 1000 Ω , dc-coupled; Inhibit impedance 650 Ω , dc-coupled; BNC connector for each.

OUTPUT

Rated output range 0.2 to 10 V positive; 12 V maximum. Output impedance ~ 2 Ω , dc-coupled, short-circuit protected; BNC connector.

ELECTRICAL AND MECHANICAL

POWER REQUIRED The Model 426 derives its power from a standard NIM bin/power supply. The power required is +24 V, 30 mA; -24 V, 49 mA; +12 V, 16 mA; -12 V, 4.9 mA.

WEIGHT

Net 0.96 kg (2.1 lb)
Shipping 1.82 kg (4.0 lb).

DIMENSIONS NIM-standard single-width module 3.43 X 22.13 cm (1.35 X 8.714 in.) per DOE/ER-0457T.

Ordering Information

To order, specify:

Model	Description
426	Linear Gate

7.5 Koinzidenzeinheit

Electronic Systems Engineering | PREP | Pool Cat...

http://cdorg.fnal.gov/ese/prep/catalog/hardware_in/

Fermilab Computing Sector

About Services Science & Computing Internal Search 

Electronic Systems Engineering

465 Triple 4-fold Coincidence

[ESE Home](#) | [Metrics](#) | [Useful Links](#) | [Vendors List](#)
[Return to LeCroy Directory](#) | [Return to Pool Catalog](#)

GENERAL DESCRIPTION

The LRS Model 465 contains three independent high-speed general-purpose coincidence units in a single-width NIM module. Each channel has four coincidence inputs and a separate veto input which accept standard negative NIM logic levels. The logic inputs may be individually enabled or disabled without altering input cabling or termination by means of front-panel pushbutton switches. With all inputs enabled, four inputs are required. Disabling the logic inputs is equivalent to reducing the number of simultaneous negative inputs required for an output. Thus, each channel may be programmed for 4-fold, 3-fold, or 2-fold logic decisions. With only one input enabled, each channel of the 465 operates as a logic fan-out.

Once triggered by signals satisfying the input coincidence requirements, the 465 generates five NIM fast logic outputs: one pair of -32 mA negative preset outputs, one -16 mA preset complementary output, and one pair of -32 mA overlap outputs. The preset outputs are continuously adjustable from less than 5 ns to greater than 500 ns by means of a front-panel multiturn potentiometer and are independent of input overlap time, amplitude, and rate. Because it is updating, it may be retriggered even before the end of an output pulse that is already present. The overlap outputs are equal in duration to the coincidence overlap and produce outputs up to the maximum input rate capability.

The front-panel fast veto input accepts standard negative NIM-level pulses. To veto the linear outputs, the veto signal must completely overlap any input coincidence; to veto the preset outputs, a prompt overlap of the leading edge of the input signal that would otherwise create the coincidence condition is required. A rear-panel bin-gate switch permits a slower (50 ns response time) inhibiting of the 465 by a clamp to ground from +4 volts through the bin gate pin of the rear power connector.

The Model 465 is packaged in a standard AEC/NIM #1 module (AEC Report TID-20893) and uses exclusively Lemo-type front panel connectors.

SPECIFICATIONS

INPUT CHARACTERISTICS

Logic Inputs:	4; Lemo-type connectors; 50 Ohm impedance; negative NIM-level input requirements; each input can be separately enabled or disabled by front-panel push-buttons.
Veto Input:	Standard negative NIM-level signal, 3.5 ns minimum width. Requires complete overlap of input coincidence for linear outputs and prompt overlap of the leading edge of the input signal that would otherwise create the coincidence condition for the preset outputs. (Veto should precede this leading edge by approximately 5 ns in this case.)
Bin Gate:	Via rear connector; clamp to ground from +4 volts inhibits; risetimes and falltimes < 50 ns.

OUTPUT CHARACTERISTICS

Preset Outputs:	3; one dual negative (quiescently 0 mA, -32 mA during output), one
-----------------	--------------------------------------------------------------------

Overlap Outputs:	positive (quiescently -16 mA, 0 mA during output). Updating. One dual negative; quiescently 0 volts, -32 mA during output; duration equal to coincidence overlap. Non-updating.
Fan-Out:	5-fold, if each output drives a 50 Ohm load.
Duration:	Continuously adjustable from less than 5 ns to greater than 500 ns by means of front-panel, screwdriver-adjustable potentiometer. Width stability: better than +/- 0.2%/°C maximum.
Output Risetimes:*	OUT: <= 2.0 ns typical (max. 2.2 ns). /OUT: <= 2.2 ns typical (max. 2.5 ns; 3.0 ns with negative output unterminated).
Output Falltimes:*	OUT: <= 2.0 ns typical (max. 2.5 ns). Slightly longer on wide output durations. /OUT: <= 2.2 ns typical (max. 2.5 ns). Slightly longer on wide output durations.

GENERAL

Logic:	2-fold, 3-fold, or 4-fold coincidences plus fan-out determined by selectively disabling logic input.
Coincidence Width:	1 ns up, determined by input pulse durations.
Rate:	120 MHz minimum.
Input-Output Delay:	13 ns preset; 8.5 ns linear.
Multiple-Pulsing:	None; one and only one output pulse of preset duration is produced each time the input conditions are satisfied regardless of the duration of the input pulses or their overlap.
Double-Pulse Resolution:	8 ns.
Packaging:	Single-width AEC/NIM module; in conformance with AEC standard; Lemo-type connectors used for all inputs and outputs.
Current Requirements:	+12 V at 65 mA +6 V at 125 mA -24 V at 5 mA -12 V at 135 mA -6 V at 640 mA


*-100 mV to -700 mV.

For assistance, contact [Help Desk](#), helpdesk@fnal.gov
Information compiled and maintained by W. Barker ; last modified on October 22nd, 2008.
(Address comments about page to ese-pcadmin@fnal.gov)

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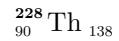2.2 Gamma Transitions and Internal Conversion Coefficients

	Energy keV	$P_{\gamma+\text{ce}} \times 100$	Multipolarity	α_K	α_L	$\alpha_M +$	α_T
$\gamma_{4,2}(\text{Ra})$	74,4 (1)	0,016 (6)	[E2]		28,9 (8)	7,89 (16)	39,5 (8)
$\gamma_{1,0}(\text{Ra})$	84,373 (3)	26,6 (14)	E2		15,9 (3)	5,8 (1)	21,7 (4)
$\gamma_{2,1}(\text{Ra})$	131,612 (4)	0,155 (8)	E1	0,195 (4)	0,041 (1)	0,013 (1)	0,249 (6)
$\gamma_{5,4}(\text{Ra})$	142,7 (1)	0,0000041 (13)	[E2]	0,280 (6)	1,396 (28)	0,50 (1)	2,18 (4)
$\gamma_{3,1}(\text{Ra})$	166,410 (4)	0,205 (16)	E2	0,225 (5)	0,704 (14)	0,256 (5)	1,185 (24)
$\gamma_{5,3}(\text{Ra})$	182,3 (1)	0,0000056 (20)	[E1]	0,090 (2)	0,0178 (3)	0,0060 (1)	0,114 (2)
$\gamma_{4,1}(\text{Ra})$	205,99 (4)	0,0201 (11)	[E1]	0,0676 (14)	0,0131 (3)	0,0042 (1)	0,0849 (17)
$\gamma_{2,0}(\text{Ra})$	215,985 (4)	0,243 (22)	E1	0,0605 (12)	0,01160 (25)	0,0038 (1)	0,0759 (15)
$\gamma_{6,3}(\text{Ra})$	228,4 (2)	0,000025 (5)	[E2]	0,125 (2)	0,182 (4)	0,065 (1)	0,372 (7)
$\gamma_{7,2}(\text{Ra})$	700,4 (1)	0,0000029 (9)	E1	0,00508 (10)	0,00084 (2)	0,000270 (5)	0,00619 (12)
$\gamma_{8,3}(\text{Ra})$	741,87 (1)	0,0000014 (4)	[E2]	0,0121 (2)	0,00330 (6)	0,00110 (2)	0,0165 (3)
$\gamma_{7,1}(\text{Ra})$	832,0 (1)	0,000014 (2)	E2+M3	0,0098 (2)	0,00240 (5)	0,00090 (2)	0,0131 (3)
$\gamma_{8,1}(\text{Ra})$	908,28 (1)	0,0000016 (5)	[M1+50%E2]	0,0203 (20)	0,0038 (4)	0,0012 (1)	0,0253 (25)
$\gamma_{8,0}(\text{Ra})$	992,65 (6)	0,0000014 (4)	[E2]	0,00720 (15)	0,00160 (3)	0,00050 (1)	0,0093 (2)

3 Atomic Data

3.1 Ra

$$\begin{aligned}\omega_K &: 0,968 \quad (4) \\ \bar{\omega}_L &: 0,452 \quad (18) \\ n_{KL} &: 0,801 \quad (5)\end{aligned}$$

3.1.1 X Radiations

	Energy keV	Relative probability
X _K		
K α_2	85,43	61,22
K α_1	88,47	100
K β_3	99,432	}
K β_1	100,13	}
K β_5''	100,738	}
		34,9
K β_2	102,89	}
K β_4	103,295	}
KO _{2,3}	103,74	}
		11,51
X _L		
L ℓ	10,622	
L α	12,196 – 12,339	
L η	13,662	
L β	14,236 – 15,447	
L γ	17,848 – 18,412	

6 Photon Emissions

6.1 X-Ray Emissions

		Energy keV	Photons per 100 disint.	
XL	(Ra)	10,622 — 18,412	8,8 (4)	
XK α_2	(Ra)	85,43	0,0172 (8)	} K α
XK α_1	(Ra)	88,47	0,0281 (12)	}
XK β_3	(Ra)	99,432	}	
XK β_1	(Ra)	100,13	0,0098 (5)	K' β_1
XK β_5''	(Ra)	100,738	}	
XK β_2	(Ra)	102,89	}	
XK β_4	(Ra)	103,295	0,00323 (16)	K' β_2
XKO _{2,3}	(Ra)	103,74	}	


6.2 Gamma Emissions

	Energy keV	Photons per 100 disint.
$\gamma_{4,2}(\text{Ra})$	74,4 (1)	0,00039 (14)
$\gamma_{1,0}(\text{Ra})$	84,373 (3)	1,17 (5)
$\gamma_{2,1}(\text{Ra})$	131,612 (4)	0,124 (6)
$\gamma_{5,4}(\text{Ra})$	142,7 (1)	0,0000013 (4)
$\gamma_{3,1}(\text{Ra})$	166,410 (4)	0,094 (7)
$\gamma_{5,3}(\text{Ra})$	182,3 (1)	0,0000050 (18)
$\gamma_{4,1}(\text{Ra})$	205,99 (4)	0,0185 (10)
$\gamma_{2,0}(\text{Ra})$	215,985 (4)	0,226 (20)
$\gamma_{6,3}(\text{Ra})$	228,4 (2)	0,000018 (4)
$\gamma_{7,2}(\text{Ra})$	700,4 (1)	0,0000029 (9)
$\gamma_{8,3}(\text{Ra})$	741,87 (1)	0,0000014 (4)
$\gamma_{7,1}(\text{Ra})$	832,0 (1)	0,000014 (2)
$\gamma_{8,1}(\text{Ra})$	908,28 (1)	0,0000016 (5)
$\gamma_{8,0}(\text{Ra})$	992,65 (6)	0,0000014 (4)

7 Main Production Modes

Th – 230(p,t)Th – 228

Th – 230($\alpha,\alpha'2n\gamma$)Th – 228

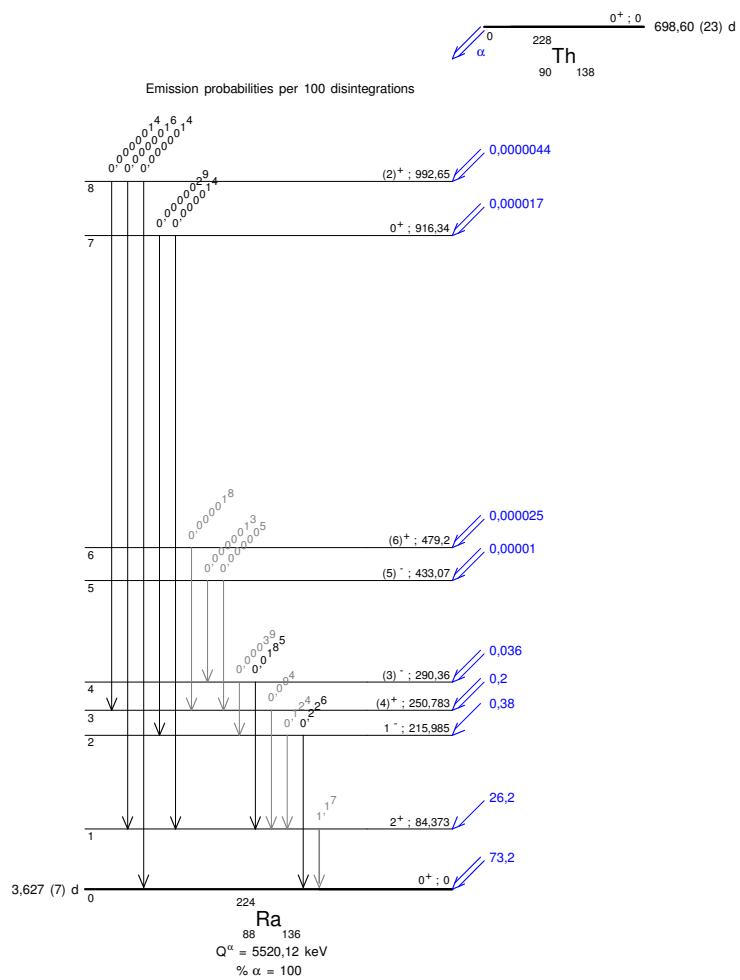8.2 Europium

1 Decay Scheme

Eu-152 disintegrates 72.1% by electron-capture and about 0.027% by emission of positrons to Sm-152 and by beta minus emission (27.9%) to Gd-152.

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
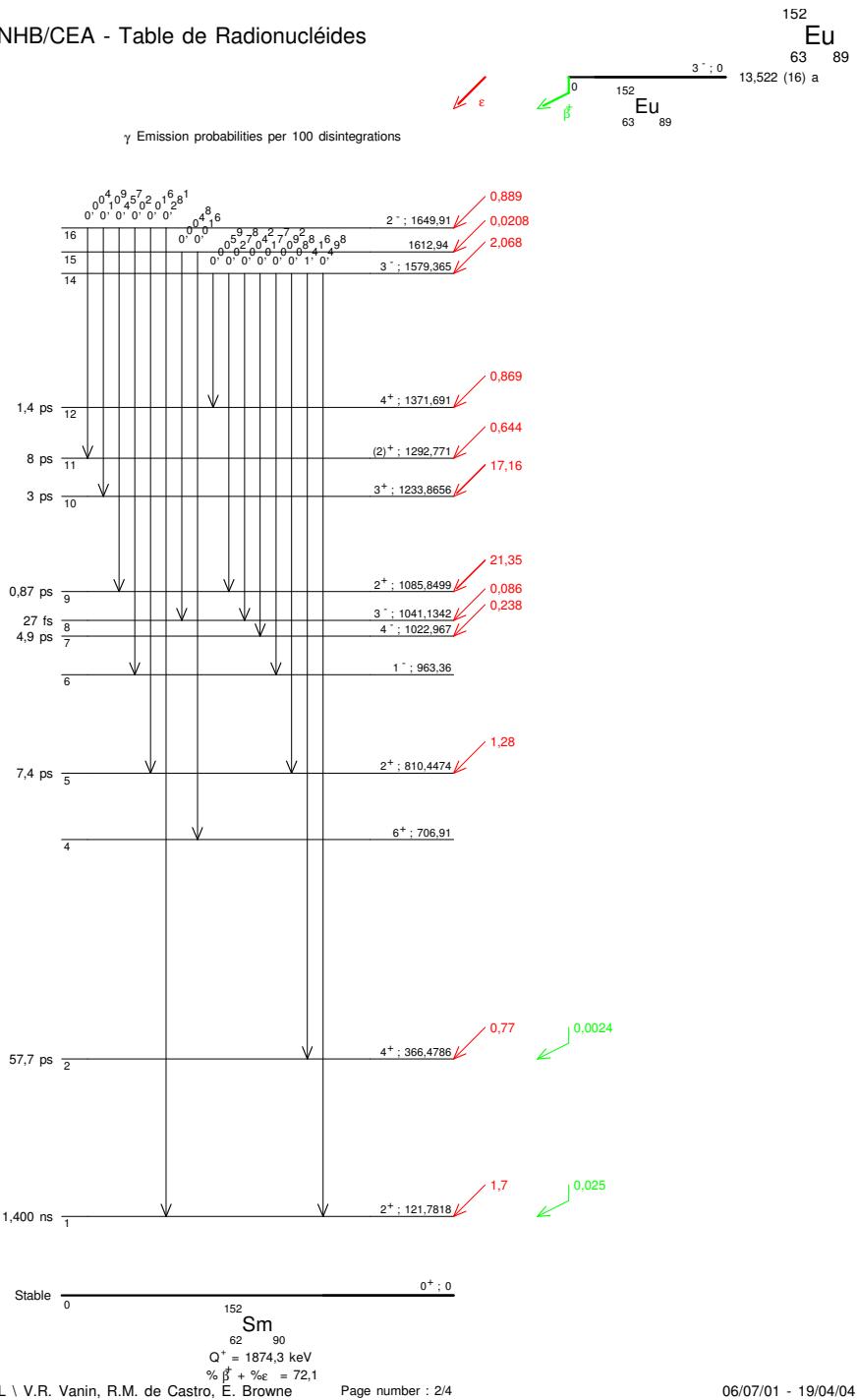
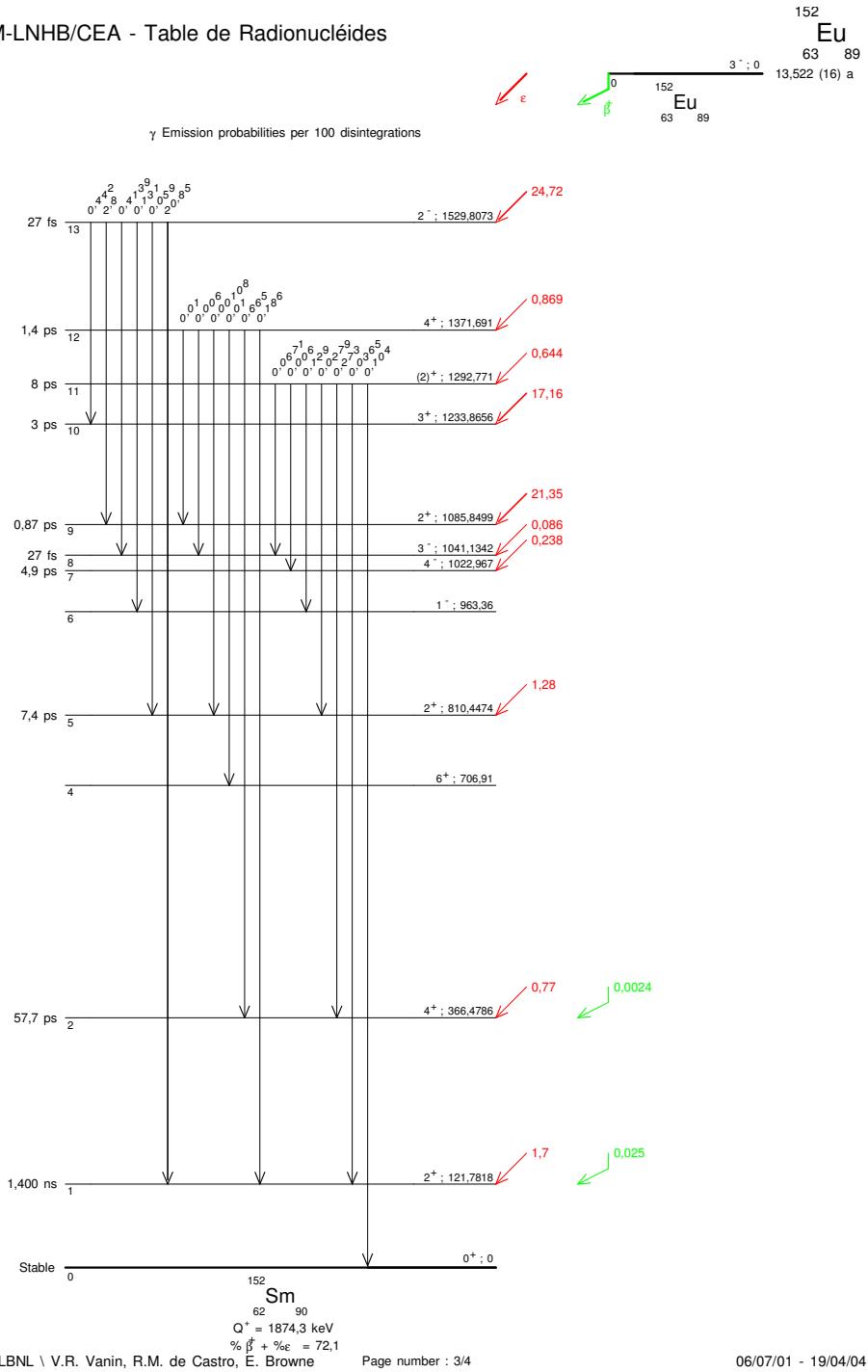
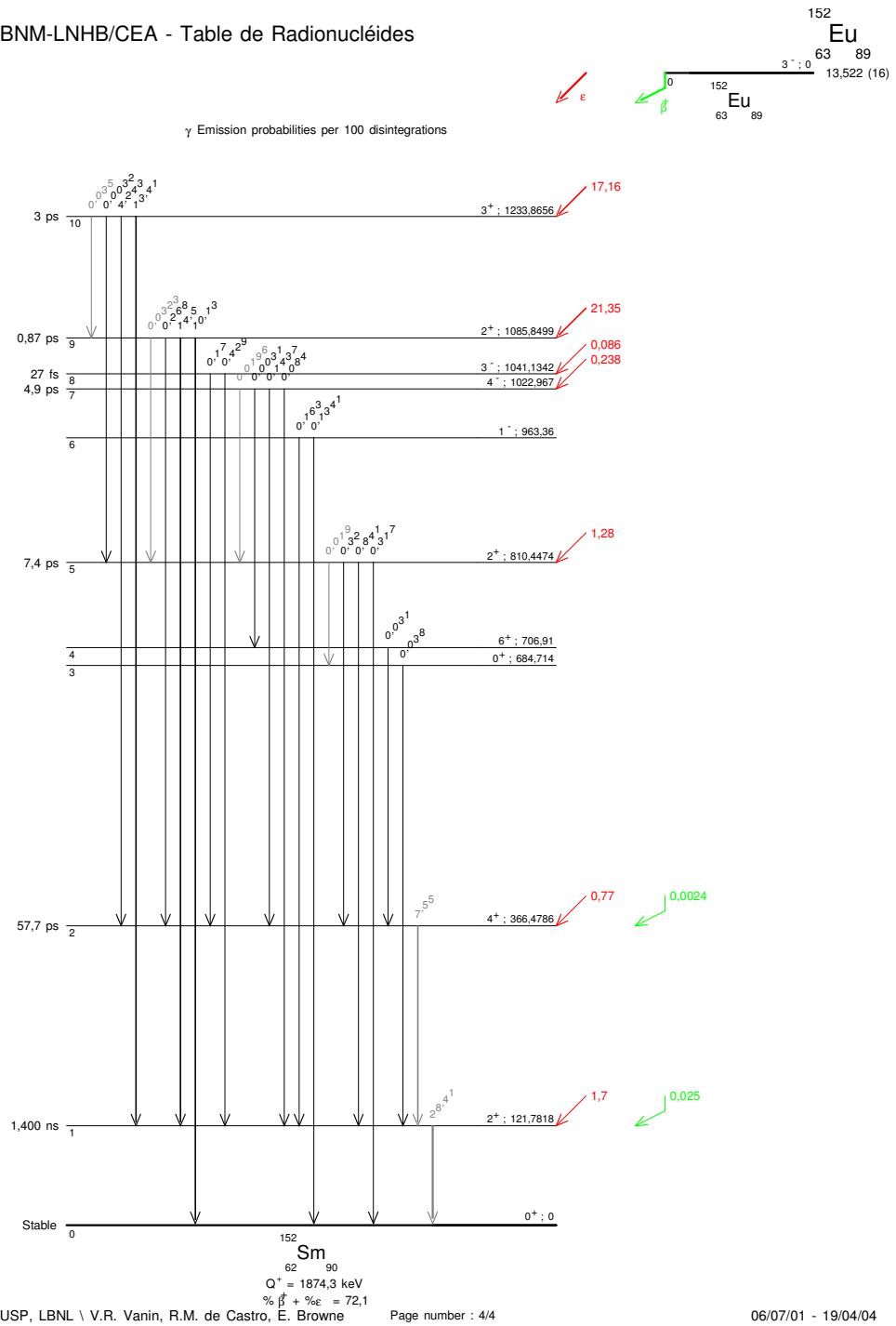
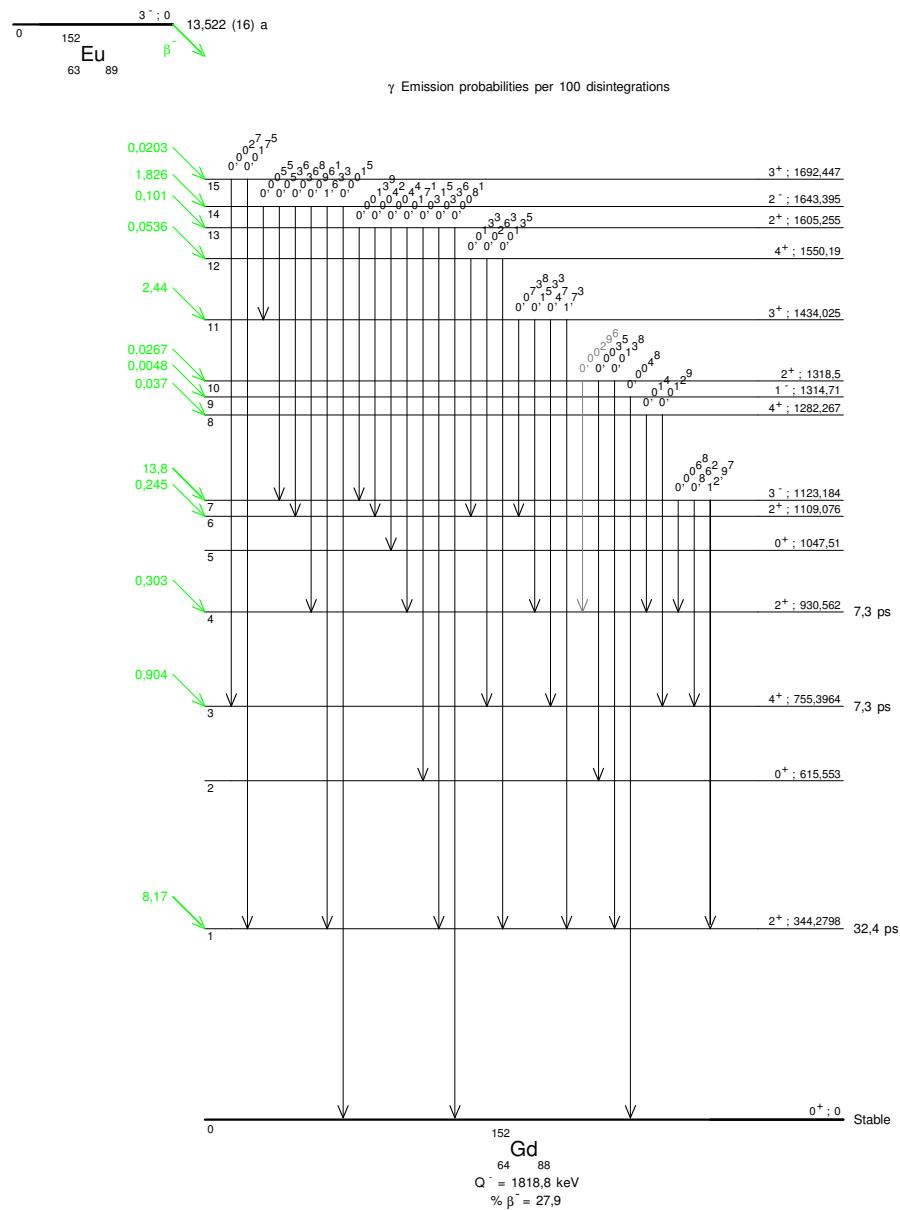
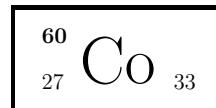
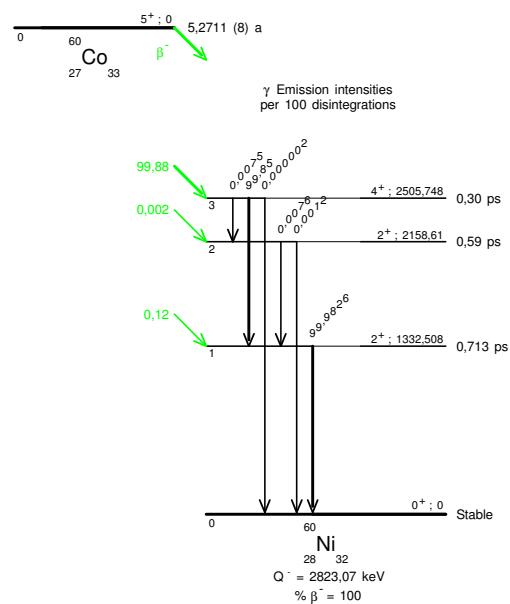
$^{22}_{11}\text{Na}$

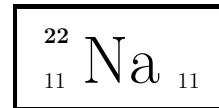1 Decay Scheme

Na-22 disintegrates predominantly to the 1275 keV level of Ne-22 . A very small fraction (0,056 %) disintegrates to the ground state of Ne-22.

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